

PRODUCT INFORMATION

Clinical Center Standardization Committee



UPDATED TUBING INSTRUCTION ABBOTT AIM PLUS AMBULATORY PUMP TUBING NOVEMBER 2003 (02)

IMPORTANT


These instructions replace

Product Information November 2003 (01)

Abbott Aim Plus Ambulatory Pump Tubing Change

Abbott Set #13732

Latex-free tubing 'Pump Set with Injector, Quick Load Cartridge, 1.2 Micron Filter, Luer Lock Adapters, Integral Anti-Valve, Nonvented, 72 inch'.

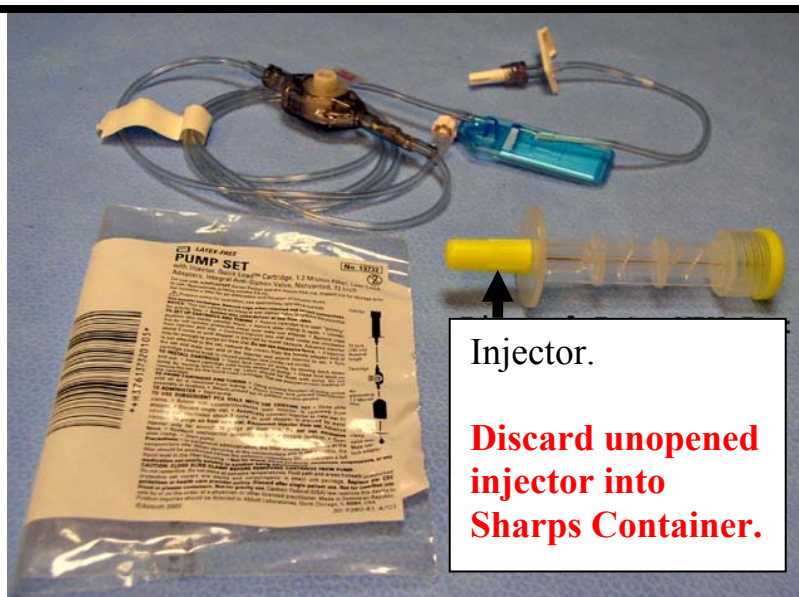

Item #: 00625 **Unit of Issue:** Each

How to find in VSC: Cartridge, Abbott,
Cartridge Infusion Provid, 00625

Instructions for Use

- ➔ Discard unopened injector into sharps container.
- ➔ The dead space is 6 ml (5.6ml priming volume). Set the pump to indicate a fill volume of 44ml for a 50 ml bag (94ml for a 100ml bag) or the patient won't get anything near the end.
- ➔ There is a disclaimer in fine print on the bag that says it may "underdeliver at low flow rates". Be careful to monitor pump.

Questions: Contact Mike Kolf, Pharmacy, email or x66551.