

# PRODUCT UPDATE

Clinical Center Standardization Committee



## UPGRADE: Kendall Kangaroo 1000ml Enteral Feeding Set December 2003 (2)

**The 1000ml enteral feeding set is now DEHP-free and has an anti-flow valve.**


Kendall Kangaroo 1000ml AFF Pump Set (Enteral Feeding Set)  
Item ID#: 02485  
Unit of Issue: Each

**How to find in VSC: Set, Kangaroo Pump 1000ml, DEHP-Free, Anti-Flow, 02485**



### Instructions:

1. Fill bag with desired amount of formula.
2. hang bag at least 18" above top of pump.
3. Add ice to pouch, if desired.
4. If using the 224/324 pump, perform the following:
  - Insert drip chamber into drip chamber guides.
  - Fill entire line with fluid by gently pulling on the black ring in a downward motion.
  - Grasp silicone tubing and gently stretch around rollers.
  - Insert retainer into retainer lock.
5. If using the PET Pump or CoNTROL Pump, perform the following:
  - Open loading arm of pump by unlatching it and rotating clockwise.
  - Insert bottom of drip chamber into lower drip chamber guide.
  - Insert top of drip chamber into upper drip chamber guide.
  - PET Pump Only—Fill entire line with fluid by gently pulling on the black ring in a downward motion.
  - Insert tubing into loading arm.
  - Close loading arm by rotating counterwise until it latches shut.

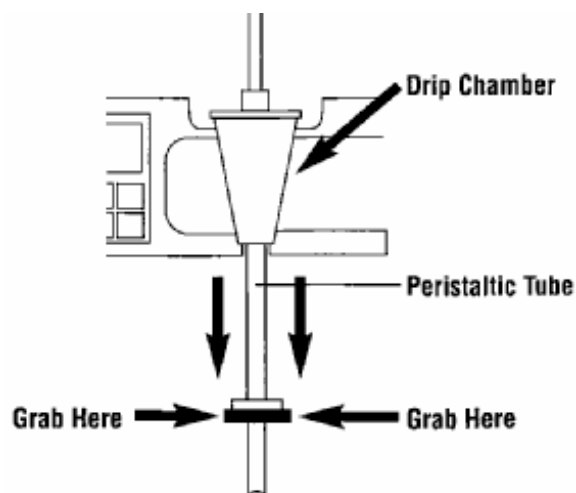

**Pull downward to Prime. Release to stop flow.**

6. Attach the pump set connector to the feeding tube.
7. Proceed with feeding.

\*\*\*\*\*