

PRODUCT UPDATE

Clinical Center Standardization Committee



April 2003 (3)

AVAILABLE FROM CHS

Blunt Filter Needle

This blunt filter needle is 18 gauge X 1 1/2". The blunt tip is a safety feature to reduce the risk for needlesticks. This needle replaces our sharp tip filter needle.

Unit of Issue: Each


Item ID#: 02754

How to find in MIS:

N—Needle—Filter

How to find in VSC:

Needle, Filter, 18G, 1-1/2

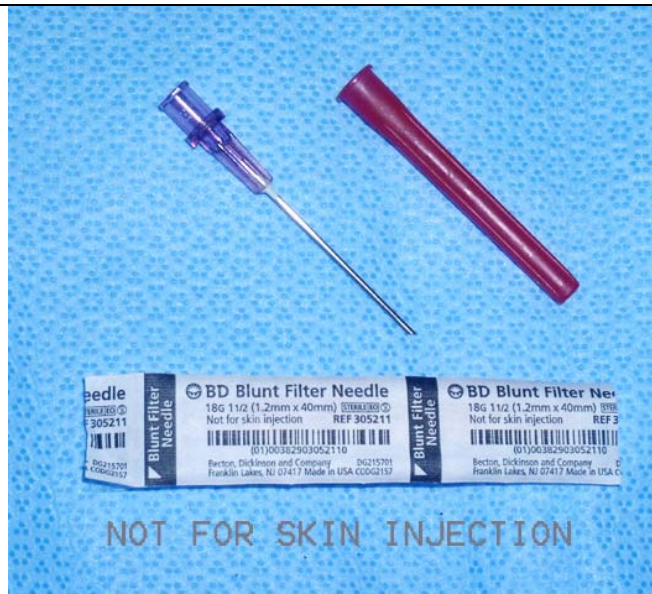

IMPORTANT: NOT FOR SKIN INJECTION.

Safety Goggle

These goggles protect the eye from aerosolized Ribavirin. The fit should be snug like swimming goggles but should not result in 'raccoon eyes'.

Unit of Issue: Each

Item ID#: 02748

How to find in MIS:

G—Goggle—Safety--Eyewear

How to find in VSC:

Goggle, Goggles, Safety, Eyewear



IMPORTANT:

- Use only during aerosolized Ribavirin therapy.
- Respiratory Therapy staff are available for consult about using the goggles.
- The goggles are assigned per person, and may be reused unless they are damaged. Clean between each use with an alcohol swab or Asepticare TB Plus II. Practice Universal Precautions at all times.