

# PRODUCT UPDATE

Clinical Center Standardization Committee



February 2003

(4)

## AVAILABLE FROM CHS

**Hemovac Compact Evacuator Kit, Middle Perforated, 10 FR, 1/8", 5" Trocar**  
Silicone wound drainage kit. 400 ml. Size = Medium. Single use.

**Unit of Issue:** Each **Item ID#: 02740**

**How to find in MIS:** D—Drain, Hemovac

**Hemovac Compact Evacuator Kit, Middle Perforated, 15 FR, 3/16", 5" Trocar**  
Silicone wound drainage kit. 400 ml. Size = Large. Single use.

**Unit of Issue:** Each **Item ID#: 02738**

**How to find in MIS:** D—Drain, Hemovac

### **Hemovac Compact Evacuator Component Kit**

Kits include 400 ml Evacuator, Evacuator Tubing, Univ. Y-Conn., 4 Silicone Adaptors.. 400 ml. Single use.

**Unit of Issue:** Each **Item ID#: 02739**

**How to find in MIS:** D—Drain, Hemovac

"For use on nursing unit when replacement evacuator is needed."

## ADDITIONAL HEMOVACS AVAILABLE FROM CHS\*

**Hemovac Compact Evacuator Kit, Middle Perforated, 19 FR, 1/4", 5" Trocar**  
Silicone wound drainage kit. 400 ml. Size = Extra-large. Single use.

**Unit of Issue:** Each **Item ID#: 01281**

**How to find in MIS:** D—Drain, Hemovac

### **Hemovac Compact Evacuator Craniotomy Kit**

Silicone wound drainage kit. 400 ml. Size = Small. Single use.

**Unit of Issue:** Each **Item ID#: 02738**


**How to find in MIS:** D—Drain, Hemovac

\* The Surgivac Hemovac has been removed from inventory.

**HEMOVAC INSTRUCTIONS: AFTER THE WOUND DRAINS ARE PLACED AND INSERTED INTO THE UNIVERSAL Y-CONNECTOR, DO THIS...**



Fully insert the evacuator tubing into the Suction Port.

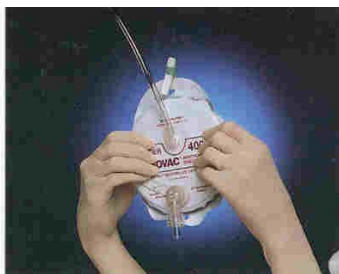

Compress the evacuator completely, making sure all the air is expelled through the Attachment Spout. Rotate the evacuator to a level position and release pressure.