


PRODUCT INFORMATION

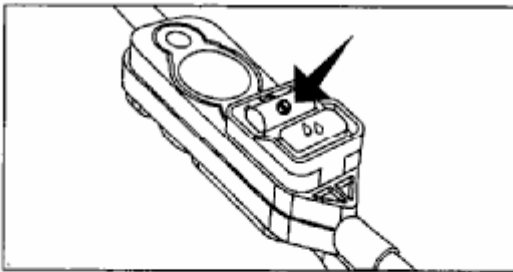
Clinical Center Standardization Committee




May 2003 (1)

GEMSTAR[®] AMBULATORY INFUSION PUMP PRIMING PUMP SETS WITHOUT A DRIP CHAMBER USING THE PURGE KEY*

- 1 Push down on the flow stop  to close the cassette.




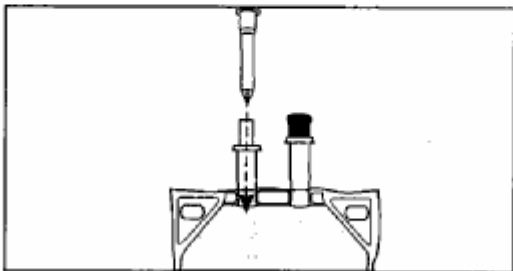- 2 Hold the container with the port toward the ceiling and insert the piercing pin.




- 3 Load the cassette.

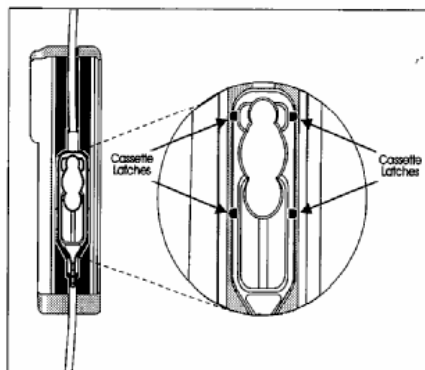CAUTION: Make sure all four cassette latches are clearly visible after the cassette is installed.



- 4 Program desired therapy.

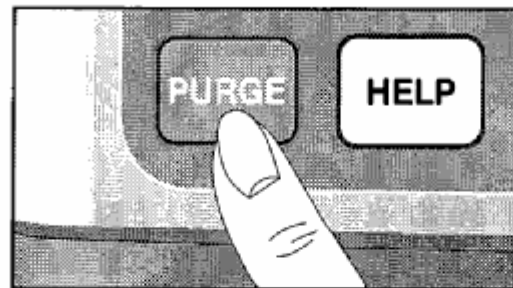- 5 Press **PURGE**.



- 6 The pump displays PRIME THE SET?

Press .

- 7 While holding the container with the port toward the ceiling, press and hold **PURGE**.

Note: Purge rate is 250 mL/hr. The pump purges for up to 2 minutes each time you press and hold

PURGE.

Continue using the **PURGE** key until fluid fills the tubing. Make sure all air is removed from the container, cassette, and tubing

Continue using the **PURGE** key until fluid fills the tubing. Make sure all air is removed from the container, cassette, and tubing

CAUTION: When priming is complete, make sure no fluid flows at the distal end of the pump set. If flow is observed, DO NOT use the pump set.

- 8 Connect the pump set to the patient-access device.

** Note: The PURGE function is only available if there is a program in the pump.
Warning: Failure to use the anti-siphon valve may result in unrestricted flow.*


PRODUCT INFORMATION

Clinical Center Standardization Committee



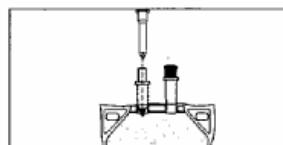
May 2003 (1)


GEMSTAR® AMBULATORY INFUSION PUMP PRIMING SETS WITH EXTENSIONS*


- 1** Push down on the flow stop  to close the cassette.




- 2** Hold the container with the port toward the ceiling and insert the piercing pin.



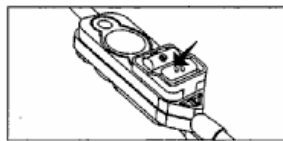- 3** Push down on the flow stop  to open (priming position) the cassette.



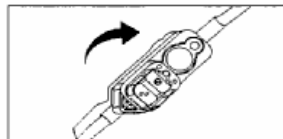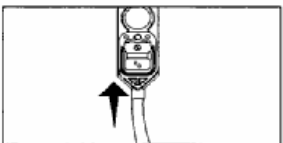
- 4** Invert the cassette approximately 45°.




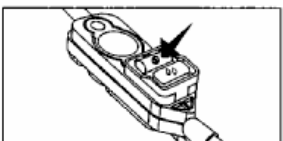
- 5** While holding the rotated cassette, squeeze or roll the container to remove air from the container and cassette.



- 6** Return the cassette to the upright position. Continue priming until fluid fills the tubing.



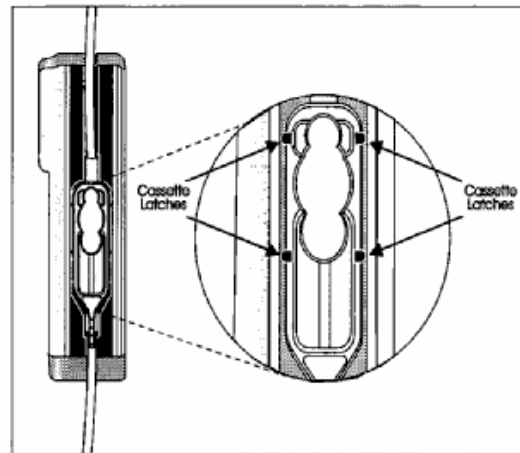
- 7** Push down on the flow stop  to close the cassette.



- 8** Connect the extension set to the pump set.

- 9** Load the cassette.

CAUTION: Make sure all four cassette latches are clearly visible after the cassette is installed.



- 10** Program the desired therapy.

- 11** Press **PURGE**.

- 12** The pump displays **PRIME THE SET?**

Press .



- 13** Press and hold **PURGE**.

Note: Purge rate is 250 mL/hr. The pump purges for up to 2 minutes each time you press and hold **PURGE**.

Continue using the **PURGE** key until fluid fills the tubing. Make sure all air is removed from the container, cassette, and tubing.

CAUTION: When priming is complete, make sure no fluid flows at the distal end of the pump set. If flow is observed, DO NOT use the pump set.

- 14** Connect the pump set to the patient-access device.

* **Warning:** Failure to use the anti-siphon valve may result in unrestricted flow.